# **PEPLYN PLUS Filters**

- liquid filters
- polypropylene



PEPLYN PLUS liquid filter cartridges are utilized for the clarification and prefiltration of a wide range of products in the fermentation and fine chemical industries.

The all polypropylene construction ensures a broad range of chemical compatibility making PEPLYN PLUS cartridges particularly suitable for the filtration of aggressive and viscous chemicals and solvents. They do not suffer from hydrolysis in aggressive solutions which would result in the contamination of the process fluid.

Extensive research has resulted in filter media with continuously graded fibre density giving progressively finer particulate retention through the depth of the media. This combined with optimized media pleating density gives PEPLYN PLUS cartridges exceptional lifetime performance.


#### Features and Benefits

- Micron rating range from 0.6 to 100 micron
- Pleated media for high flow rates and long life
- Graded density for excellent particle retention
- Wide range of end caps to provide retrofitting of existing systems
- All polypropylene construction




Note: PEPLYN is a registered trademark of Parker Hannifin Corporation.

## **Performance Characteristics**

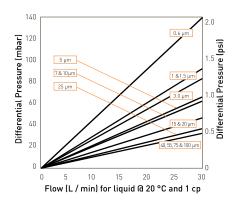

For K size for a given flow rate multiply  $10^{\circ\circ}$  size differential pressure by 2

10" Size (250 mm) Cartridge

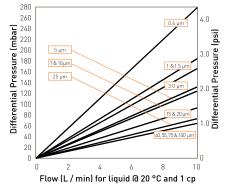

For A size for a given flow rate divide B size differential pressure by 2 For E size for a given flow rate multiply B size differential pressure by 2

B Size (65 mm) Cartridge

# **Specifications**

#### Materials of Construction

Filtration Media: Polypropylene
Upstream Support: Polypropylene
Downstream Support: Polypropylene
Inner Support Core: Polypropylene
Outer Protection Cage: Polypropylene
End Caps: Polypropylene

End Cap Insert (if applicable): 316L Stainless Steel\*
\*Not available in B & L endcap variants

■ Standard o-rings/gaskets: Silicone / EPDM

#### **Biological Safety**

Materials conform to the relevant requirements of 21CFR Part 177 and current USP Plastics Class VI - 121 °C and ISO10993 equivalents.

#### **Recommended Operating Conditions**

Up to 70 °C (158 °F) continuous operating temperature and higher short-term temperatures during CIP to the following limits:

Temperature °C °F		Max. Forward dP (bar) (psi)		
-				
20	68	5.0	72.5	
40	104	4.0	58.0	
60	140	3.0	43.5	
80	176	2.0	29.0	
90	194	1.0	14.5	
>100 (steam)	>212 (steam)	0.3	4.0	

## Effective Filtration Area (EFA)

10" (250 mm) Up to 0.79 m<sup>2</sup> (8.50 ft<sup>2</sup>)

#### Cleaning and Sterilization

PEPLYN PLUS cartridges can be repeatedly steam sterilized in situ or autoclaved at up to 135 °C (275 °F). They can be sanitized with hot water at up to 90 °C (194 °F) and are compatible with a wide range of chemicals.

#### **Retention Characteristics**

The retention characteristics of PEPLYN PLUS have been determined by a single-pass technique using suspensions of ISO 12103 Part 1 A2 Fine and A4 Coarse test dust in water.

Media Code	Micro >99.99% 10000	n Rating a 99.98% 5000	t Various E 99.90% 1000	fficiencies 99% 100	90% 10
.60	0.60	0.57	0.54	0.32	0.20
1.0	1.00	0.95	0.90	0.70	0.50
1.5	1.50	1.40	1.10	0.80	0.60
003	3.00	2.80	1.80	1.00	0.70
005	5.00	4.70	4.50	3.50	1.00
007	7.00	6.70	6.30	4.50	2.50
010	10.00	8.00	7.00	4.80	2.80
015	15.00	12.00	10.00	7.20	4.50
020	20.00	16.00	14.00	10.00	6.00
025	25.00	20.00	17.00	12.00	7.00

## **Ordering Information**

## Cartridges

